



Guided Surgery



The adinguide™ was created with you in mind, bringing the guided surgery world to your clinic.

The kits were specifically designed to support effortless procedures, accurate and predictable implant placement.

The ActiveFlow™ Irrigation Technology is a unique design that forces coolant through the guide to ensure irrigation reaches the bone and prevents overheating. The adinguide™ drills are designed with carved channels which allow saline coolant from your hand-piece to stream through the sleeve directly into the osteotomy site while drilling.

Our uniquely designed drill kit is intuitive and easy to use by implant diameter and length.

The Tissue Punch and Pilot Drill allow for easy flapless access.

Keyless procedure enables easy handling and full visibility of surgical field.

Self-centering drills – controlled by a 5.5mm diameter cylinder at the top of each drill centers the drill in the guide's sleeve.

Self-centering cylinder engages the sleeve before the drill engages the osteotomy in order to ensure correct positioning of the drill.

Built-in stoppers – ensure precise and accurate drilling to the desired depth.



Guided Surgery Kit

Implants Drills | GSK001



For full guided procedure:



RS Touareg™ S/OS
Implant Platform



GS Tools Kit RS
GSK003



5.5mmD Sleeve
GS0055



GS5055 Tissue Punch 5.0mmD, S-5.5



GS0625 Pilot Drill S-5.5, 2.5/1.8mmDx6.0mml



GS0825 Drill S-5.5, 2.5/1.8mmDx8.0mml



GS1025 Drill S-5.5, 2.5/1.8mmDx10.0mml




GS1125 Drill S-5.5, 2.5/1.8mmDx11.5mml



GS1325 Drill S-5.5, 2.5/1.8mmDx13.0mml



GS1625 Drill S-5.5, 2.5/1.8mmDx16.0mml

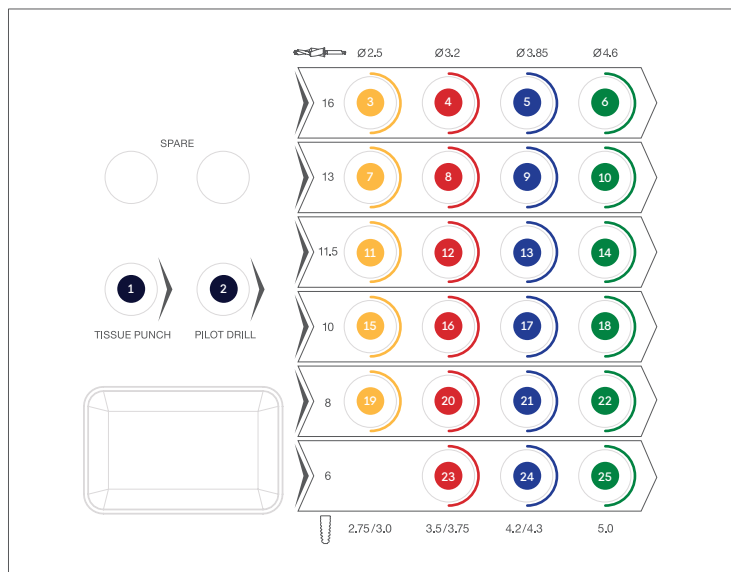

GS0632 Drill S-5.5, 3.2/2.5mmDx6.25mml



GS0832 Drill S-5.5, 3.2/2.5mmDx8.0mml



GS1032 Drill S-5.5, 3.2/2.5mmDx10.0mml



GS1132 Drill S-5.5, 3.2/2.5mmDx11.5mml



GS1332 Drill S-5.5, 3.2/2.5mmDx13.0mml



GS1632 Drill S-5.5, 3.2/2.5mmDx16.0mml



GS0638 Drill S-5.5, 3.85/3.2mmDx6.25mml



GS0838 Drill S-5.5, 3.85/3.2mmDx8.0mml



GS1038 Drill S-5.5, 3.85/3.2mmDx10.0mml



GS1138 Drill S-5.5, 3.85/3.2mmDx11.5mml



GS1338 Drill S-5.5, 3.85/3.2mmDx13.0mml



GS1638 Drill S-5.5, 3.85/3.2mmDx16.0mml



GS0646 Drill S-5.5, 4.6/3.85mmDx6.25mml



GS0846 Drill S-5.5, 4.6/3.85mmDx8.0mml



GS1046 Drill S-5.5, 4.6/3.85mmDx10.0mml



GS1146 Drill S-5.5, 4.6/3.85mmDx11.5mml



GS1346 Drill S-5.5, 4.6/3.85mmDx13.0mml



GS1646 Drill S-5.5, 4.6/3.85mmDx16.0mml

Guided Surgery Kit

RS Touareg™ Implants Tools

GSK003



For full guided procedure:




RS Touareg™ S/OS
Implant Platform



GS Drill Kit
GSK001



5.5mmD Sleeve
GS0055

1

GS0032
RS Handpiece Fixture Mount



4

GS1917
Anchoring Drill S-2.0, 2.0/1.7mmDx19mmL



2

GS0034
RS Grip Fixture Mount



5

GS1915
Anchoring Drill S-2.0, 2.0/1.5mmDx19mmL




3

GS0029
RS Screw Fixture Mount (qty 4)



6

GS0028
Handpiece Guide Fixation Screw Driver

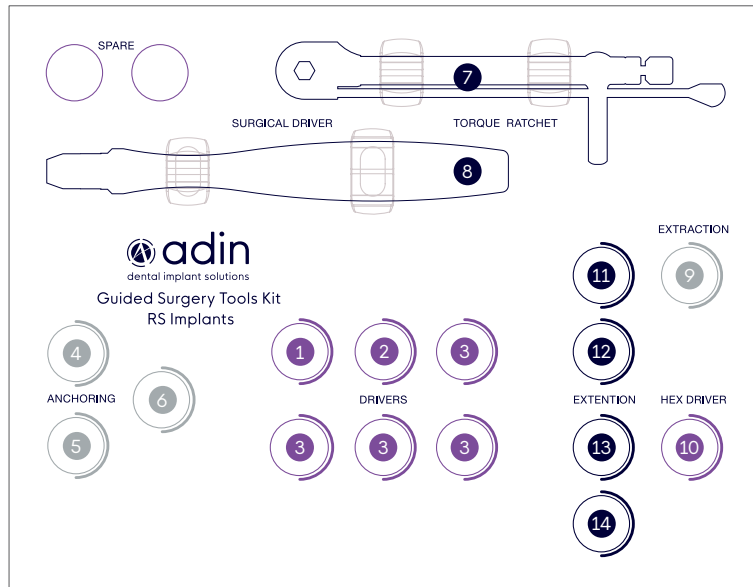

7

GS6111
Surgical Torque Ratchet 35-100 Ncm for hex 4mm



10

RS6094
1.27 Hex Hand Driver-Medium



13

GS0021
Handpiece Socket L10.5



8

GS0027
Guide Fixation Screw Driver



11

GS0023
Manual Socket L10.5



14

RS9027
Drill Extension



9

GS0035
RS Extraction Tool



12

GS0024
Manual Socket L14.5



Additional Products

Sold separately

Accessories



2.0mmD Sleeve

Catalog No. GS0020



5.5mmD Sleeve

Catalog No. GS0055



Guided Vertical Fixation Screw S-3.3

Catalog No. GS0006



Guided Vertical Fixation Screw S-5.5

Catalog No. GS0007



2.0D 19mmL Anchorage Screw

Catalog No. GS1920

| Purchase Options

Catalog No.	Description	Package contents	
GS0006	Guided Vertical Fixation Screw S-3.3		x 3 pack
GS0007	Guided Vertical Fixation Screw S-5.5		x 3 pack
GS1920	2.0mmD Anchoring Screw		x 3 pack
GS0055	5.5mmD Sleeve		x 1 pack
GS0055x3			x 3 pack
GS0055x10			x 10 pack
GS0020	2.0mmD Anchoring Sleeve		x 3 pack



18mm Drills



Drill S-5.5, 2.5/1.8mmDx18.0mmL

Catalog No. GS1825



Drill S-5.5, 3.2/2.5mmDx18.0mmL

Catalog No. GS1832



Drill S-5.5, 3.85/3.2mmDx18.0mmL

Catalog No. GS1838



Drill S-5.5, 4.6/3.85mmDx18.0mmL

Catalog No. GS1846



Drill S-3.3, 2.5/1.8mmDx18.0mmL

Catalog No. GSN1825



Illustrations are not to scale
 Availability of the products in some places is limited and depends on local regulations.
 For more information, please contact your local representative.